

1 litre makes up to 125 litres of
disinfectant or 250 litres of cleaner

New PULI-JET PLUS

Disinfectant for Dental Vacuum and Amalgam Separators




- ① Unscrew red cap slightly.
- ② Squeeze bottle to distribute liquid to top container.

- ③ To mix, fill the 'Pulse Cleaner Bucket' with warm water (**approx temp 50°C**) up to the third marker, then pour the Puli-Jet Plus into the 'Pulse Cleaner bucket'.

General Cleaning: We recommend that you cleanse the aspiration system **every evening**, at the **end of work**. For this mix 10ml of Puli-Jet Plus with 2.5 litres of warm water.

Disinfecting: We recommend to disinfect the aspiration system **once a week** with a higher concentrate of Puli-Jet Plus. For this mix 20ml of Puli-Jet with 2.5 litres of warm water.

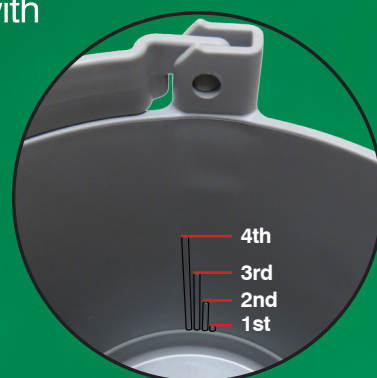

- **Non-foaming** concentrate
- Without **aldehydes** and **phenol**
- For disinfecting and simultaneous deodorization, cleaning and maintenance of dental aspirators and amalgam separators (wet, semi-wet & dry systems)
- Please always refer to the manufacturers instructions when using Puli-Jet Plus

**LASTS
LONGER
THAN MOST
OTHER
BRANDS**